



GEMaC


Groupe d'Étude de la Matière Condensée

TRANSPORT MEASUREMENTS

Aligned four points probe for resistivity measurements

Non-destructive method to measure resistance (resistivity) of samples at room temperature.

Two sets: 0.5 mm and 2 mm distance between probes (WC tips), with XY stages

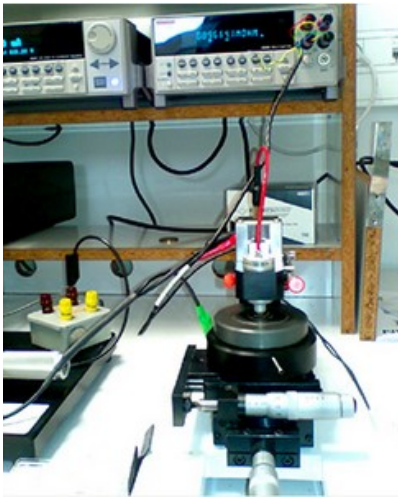

I-V current-tension and C-V capacitance measurements

A Keithley 4200 “Semiconductor characterization system” equipped with 5 probes allows:

- » I-V characterization. Measures 2 and 4 points probe resistivity at room temperature. From femtoA to A currents, from nanoV to 200 V
- » Capacitance – Voltage measurements. From attoFarads (aF) to nanoFarads (nF) range at frequencies from 10 KHz to 10 MHz. With two Kelvin probes
- » In general, any tests of semiconducting IC with an “Labview-type” commercial program



Measure capacitance from attofarads (aF) to nanofarads (nF)

Test at frequencies from 10kHz to 10MHz

Resistivity and Hall effect measurements

We have developed two “Resistivity and Hall effect measurements” set-ups in Van der Pauw configuration.

The first set-up is dedicated to measurements at room temperature of middle range impedance ($1e+8$ Ohms) samples up to 0.8 T magnetic field.



The second “high impedance-high temperature” ($1e+10$ Ohm) set-up is equipped with a 1.6 T electro-magnet in 80 K - 880 K temperature range.



We study various semiconductors and oxides from their highly conductive up to semi-insulating state.